


ELDONCARD®2511 for ABO & RhD forward blood typing

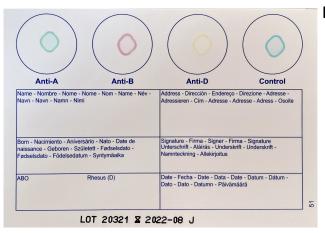

Recommended for:

- Medical Market Blood Banks, Hospitals, Clinics etc.
- ☑ Emergency & Rescue NGO's
- ☑ Donor Recruitment
- ☑ Rhesus tests for Pregnant Women
- ✓ Educational Market
- ✓ Health / Home Market
- ☐ Bed side test

In the circles, dried blood typing reagents are found. The card is used for the detection of the blood type of one person (A, B, AB, or O and D (RhD positive or negative) in each of the 4 circles of one drop of water and one drop of blood from the testperson. Control field to detect false possible false positives. The procedure is described in detail in the instructions for use.

Available as:



ELDON BAG Contains:

 25 ELDONCARD 2511 25 Blood Group Determinations



- 100 ELDONCARDS 2511 100 Blood Group Determinations
- ELDON BIG BAG Contains:*
 - 50 or 100 ELDONCARDS 2511 50/100 Blood Group Determinations

Costume made Kits for special use NGO/Doctors/

Military

*Each EldonCard is in an aluminum envelope.

Labelling: CE & FDA approved and a large number of other

countries

Storage: 5-37 Celsius / 41-99 Fahrenheit – NO COOLCHAIN

Expiration: 24 month from first day of production

Accessories:

In general, accessories are available in numbers matching the number of Blood Group Determinations.



Safety Lancet



Steel Lancet



EldonStick



Alcohol Swabs



EldonFoil



Pippette

General: To perform one Blood Group Determination 2511 you need 1 lancet, 4 EldonSticks, 1 EldonFoil, 1 pipette and 1 alcohol swab.

Contacts: Kasper Juul Hedegaard - ki@eldoncard.com - +4542138115

